

PRODUCT CODE:

MS04010004



velo

Serial Number

INSTRUCTION MANUAL

CONTENTS

1. INTRODUCTION	2
2. IMPORTANT PRECAUTIONS	3
2.1 BEFORE DRIVING SCOOTER	
2.2 CAUTIONS WHEN DRIVING	
2.3 CAUTIONS WHEN FOLDING/UNFOLDING SCOOTER	
2.4 CAUTIONS WHEN DRIVING ON INCLINE	
2.5 CAUTIONS WHEN CLIMBING	
2.6 ELECTROMAGNETIC INTERFERENCE AND WARNINGS	
2.7 SAFETY WARNING AND INSTRUCTION LABELS	
3. IDENTIFICATION OF PARTS	12
4. INTRODUCTION & OPERATION OF SCOOTER	13
4.1 IDENTIFICATION & OPERATION OF CONTROL PANEL	
4.2 OPERATION OF MAIN SWITCH	
4.3 OPERATION OF THROTTLE	
4.4 OPERATION OF FREE-WHEELING	
4.5 TILLER HEIGHT ADJUSTMENT	
4.6 SEAT ADJUSTMENT / DISMANTLE	
5. FOLDING / UNFOLDING YOUR SCOOTER	20
6. BATTERY GAUGE & CHARGING THE BATTERIES	23
6.1 DISASSEMBLING BATTERY BOX	
6.2 CAUTIONS FOR BATTERY & CHARGER	
6.3 BATTERY GAUGE	
6.4 OPERATING THE CHARGER	
7. CARE AND MAINTENANCE	29
7.1 INSPECTION CHECKS	
7.2 CLEANING AND MAINTENANCE	
8. TROUBLE SHOOTING	32
9. OTHER INFORMATION	34
9.1 OPTIONS	
10. TECHNICAL SPECIFICATIONS	36
11. WARRANTY	37

1. INTRODUCTI

Thank you and congratulations on purchasing your new Li-Tech Velo mobility scooter.

It is designed to provide the transportation ability indoors and outdoors for person whose ability to walk is impaired, but who are still in terms of their eyesight and physically and mentally able to operate an electric scooter.

We pride ourselves on providing safe and comfortable products. Our goal is to ensure your complete satisfaction. We sincerely hope you enjoy your Li-Tech Velo mobility scooter.

Please read and observe all warnings and instructions provided in owner's manual before your operate the various functions of this scooter. Also, please retain this booklet for future reference.

If you have any question, please contact:



CareCo (UK) Limited
1 Turing Court
Great Notley
Braintree
CM77 7AT
England

Tel: 0333 015 5000
www.CareCo.co.uk

INFORMATION OF EUROPEAN REPRESENTATION



MedNet EC-REP GmbH
Borkstrasse 10, 48163 Muenster,
Germany



MediMap Ltd, 2 The Drift,
Thurston, Suffolk, IP31 3RT

2. IMPORTANT PRECAUTIONS

2.1 BEFORE DRIVING SCOOTER

- Before take the first trip with mobility scooter, you should familiarize yourself well with the operation of mobility scooter and with the operation elements. Please take your time to read this introduction booklet.
- Before driving, please evaluate the personal condition, and fully understand the operation of mobility scooter.
- Before charging the scooter, please release the Circuit Breaker on battery box first to activate the system (see section 3)
- Scooter is shipped with 30% of battery charge only, batteries should be fully charged before your journey for maximum performance and longevity (see section 6).
- Scooter is equipped with certain type of battery (lithium or optional lead-acid), if you want to change it, please contact dealer to reset the battery system such as charging port and controller setting. Only one type of batteries can be used on scooter, do not mix the batteries and charge, or it may result dangers
- The maximum load of the scooter is 115KG. Do not exceed the maximum permissible load. Exceeding the max weight rating may result in injury to yourself.
- The mobility device is only designed for use by a single occupant whose maximum weight does not exceed the maximum permissible load of the device. Never use the mobility device to transport more than one person (including children).
- You should not assemble, maintain, and operate the mobility scooter before you read this instruction booklet.
- Do not attempt to modify or dismantle the battery box or carry out any maintenance work that is not described in this user manual.
- Do not change, modify, remove any parts (seat, shroud) from products especially safety protected parts such as anti-tippers.
- Completely deflate the tyres (air tyres) before dismantling the rim or attempting.
- Materials and assemblies of scooter are flame resistant.
- Always keep battery box installed on scooter to avoid short-circuit, do not attempt to remove battery box at all possible.
- Always keep the tiller screw knob fasten at whole time except when adjusting the tiller height.
- Do not sit or stay on scooter when scooter is on moving transport vehicles to avoid dangers
- Turn the power (key) off before getting on or off. (see section 4-2)
- Do not drive your scooter if you are under the influence of alcohol or medication that may affect your ability.
- Do not turn off the power while driving. Doing so will bring the scooter to an abrupt stop and danger.
- Please observe all relevant rules and regulations pertaining to pedestrians and road users at all times when you are driving the scooter.

2. IMPORTANT PRECAUTIONS

- Mobility device may only be used on the traffic routes for which it is approved in accordance with the relevant national legislation.
- Always be aware of pedestrians and situations which might require extra care when using your scooter on public walk ways and footpaths.
- Avoid driving scooter in dark or in conditions of bad weather and visibility for safety



- Do not turn the power on before you get in and sit securely on the seat. Be certain that the power is turned off when get in, get out. This will eliminate the possibility of accidentally activating the wigwag controls and causing injury.
- Keep your weight toward the middle of the deck. Putting most or all of your weight on the edge of the deck may cause an unstable condition.
- Do not use hose pipe or splash water directly onto the scooter or store the scooter in a humid condition such as rain or show as this may cause damage to electronics.
- Do not remove safety protected parts such as anti-tipper from scooter.
- Do not use the unauthorized parts such as battery or charger on the scooter.
- Use the accessories supplied from original dealer only.
- Do not exceed the maximum permissible load of scooter or it may cause unbalance.
- Please stop operating the scooter and contact your dealer if you find the sign of scooter malfunction.
- Please pay more attention when driving the scooter, the emergency stop may encounter (ex. out of power, EMI interference)
- Keep the RF sensor key (Remote) carefully. Please contact dealer if the key is lost.



Practice Tips

- If you are new to driving a scooter, it is a good idea to practice in a clear, safe space on an flat level surface. (I.e. Park, Playground).
 - Basic functions to practice: Wig wag accelerate /Wig wag release /Stop / Reverse /Turn /Ramp proceeding.
 - Set the speed control to its lowest speed, slightly increase the speed when you are getting familiar with the scooter.
-

2. IMPORTANT PRECAUTIONS



- Practice operating your scooter in the presence of an attendant. Remember that only with practice will you become a competent driver.
 - Practice these basic functions until you feel that you have control of your scooter.
-



Accessories Contents for Your Scooter :

- Auto-folding version : control remote
-

2.2 CAUTIONS WHEN DRIVING

- Please do the daily check before your journey always. (See section 7-1)
- Do not stretch your body out on the scooter.
- Do not use the mobile phone and wireless mobile devices while driving.
- Do not use the batteries of mobility scooter to charge any other electric devices, except the accessories from original manufacturer.
- Do not drive on a slope which is over the limit. (See section 2-4)
- Please bypass, do not drive to across the obstacle over the limit (See section 2-5)
- Do not attempt to drive the scooter under below hazardous road condition. If necessary, please drive it in the presence of an attendant:
 - Dark, bad weather such as in rain, snow, mist, storm, and poor visibility.
 - Narrow paths, level crossing, and uneven terrain.
 - Wet grass, sand, gravel, or any other potentially hazardous condition.
- Please stop operating the scooter if the batteries have drained, continuous operation may damage the scooter.
- Do not drive when the seat lock is positioned in folding point (red). Scooter can be operated only when the seat lock is in the fixation point (green). (see section 4-6)
- When turning, reduce your speed and maintain a stable center of gravity. This greatly reduces the possibility of a tip or fall.
- Do not bounce or jam on the footplate.
- Please always set the speed control to its lowest speed when you start to drive, and observe all relevant rules and regulations pertaining to pedestrians and road users at all times.

2. IMPORTANT PRECAUTIONS



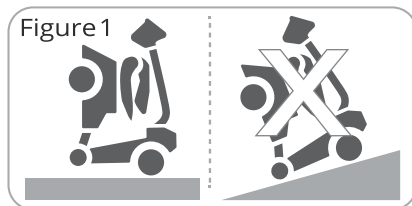
- Reduce speed when descending to prevent any danger. Increase speed when ascending (see section 2-4)
 - Do not turn off the power while driving.
 - Do not operate on staircases or escalators.
 - Do not attempt to climb curbs greater than limitation show on Technical Specification.
 - Do not make S turns or sharp turns while driving. This may cause the scooter to tip.
 - Never place the scooter in freewheel mode (Neutral) when on any sort of an incline or decline to avoid roll-away. .
 - When the scooter is parked, the lever for engaging and disengaging the motors must be locked firmly into the "DRIVE" position. (see section 4-4)
 - Do not conduct any adjustment such as tiller / light setting / control panel setting while driving.
 - Do not pull the seat springs while driving to avoid dangers. (see section 4-6)
-

2.3 CAUTIONS WHEN FOLDING/UNFOLDING SCOOTER


- When folding or unfolding, you must place the scooter on a flat ground to avoid tip-over, roll-away, or dangers. Never operating the folding/unfolding process on a slope
- Please read through the section 5 for folding / unfolding operation before you operate the scooter.
- Turn off the power before you fold / unfold the scooter. (see section 4-2)
- Folding/unfolding actions will stop if any malfunctions are detected.(see section 8)
- Please make sure the lever for engaging and disengaging the motors must be locked firmly into the "DRIVE" position before you fold / unfold the scooter. Never proceed the folding/unfolding under freewheel mode "Neutral" position (see section 4-4)
- Scooter can be charged when it's folded or unfolded. However, scooter will not respond to folding / unfolding commands when charger is plugged. To operate the folding / unfolding, please unplug the charger from the scooter (see section 6-4).
- Please remove all the accessories, such as armrests or bags out before folding / unfolding the scooter
- Please observe environment surround you for safety before folding/unfolding the scooter.
- When the folding or unfolding is processing, please pay attention to the movement of tiller for your safety.

2. IMPORTANT PRECAUTIONS

- When the folding or unfolding is processing, do not sit or interrupt the process or it may cause damage or danger to the scooter and yourself.
- Please double check and make sure the procedure of folding/unfolding is completed before you get on the scooter.
- Do not turn on the power when the scooter is folded.
- Scooter can be rolled by front wheels and supporting castors when it's folded (refer to Right picture).
- Do not sit or stay on folded/unfolded scooter when scooter is on moving transport vehicles.
- Do not store the scooter in a humid condition such as rain or show as this may cause damage to electronics. (see section 7-2)
- The brake does not apply when the scooter is in folded and standing position. Do not place the standing folded scooter on a slope to avoid roll-away and danger (Fig.1)
- When folding or unfolding, you must place and lay down the scooter on a flat ground (Fig. 2)



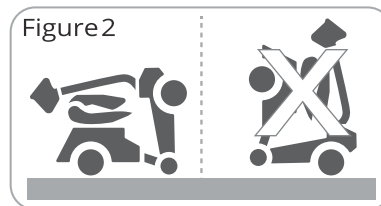Do not place scooter on a slope when scooter is folded and standing.



Place and laydown the scooter before unfolding.

2.4 CAUTIONS WHEN DRIVING ON INCLINE

- The scooter has been rated to a maximum climbable height, obstacle height, and gap. (see section 10)
- Never drive on a slope that exceeds the rated slope.
- For maximum stability, lean forward in your body while proceeding up ramps, inclines, curbs, or any low rise. (please refer to following picture)



On flat road



On slope



2. IMPORTANT PRECAUTIONS

- The weight capacity limit at different ramp degree (please refer to above picture).
- Your scooter's ability to travel up inclines is affected by your weight, your scooter's speed, your angle of approach to the incline, and your scooter setup.
- Please avoid to driving on a long ramp or any uneven terrain to prevent any danger from motor defected.



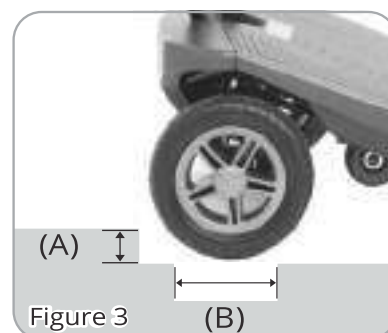
- When driving down a ramp or uneven terrain, keep the scooter's speed adjustment set to the slowest to ensure a safely controlled driving.
 - If the speed is too fast, release the throttle control lever, let the scooter stop. When you feel that you again have control of your scooter, push the throttle control lever forward and continue safely driving.
 - To prevent any danger, do not make turns on ascending, descending ramp.
 - The batteries voltage normally will go up when driving on descending road. If the battery voltage becomes too high, the over-voltage protection will be activated by slowing the speed till the scooter stops. (error light indicator will flash 3 times). Please pull over the scooter to the safe area, release the wigwags and restart the scooter again.
-

2.5 CAUTIONS WHEN CLIMBING

- The maximum height of obstacle and curb that scooter can climb is up to 5 cm (A).
- The maximum gap that scooter can drive over is 10 cm (B).
- When driving scooter on ramp, adjust body center of gravity for scooter stability (see section 2-4).



- Even though your scooter has been rated with maximum climbable, we suggest not driving through an obstacle.
-



2. IMPORTANT PRECAUTIONS



- It is very important that you read this information regarding the possible effects of Electromagnetic Interference on your mobility scooter.
-

2.6 ELECTROMAGNETIC INTERFERENCE AND WARNINGS

Mobility scooters may be susceptible to electromagnetic interference (EMI), which is interfering electromagnetic energy (EM) emitted from sources such as radio stations, TV stations, amateur radio (HAM) transmitters, two-way radios, and mobile phones. The interference (from radio wave sources) can cause the mobility scooter to release its brakes, move by itself, or move in unintended directions. It can also permanently damage the mobility scooter control system. The intensity of the interfering EM energy can be measured in volts per meter (V/m). Each mobility scooter can resist EMI up to certain intensity. This is called its "immunity level." The higher the immunity level, the greater the protection. At this time, current technology is capable of achieving at least a 20 V/m immunity level, which would provide useful protection from the more common sources of radiated EMI. The immunity level of this mobility scooter model is not known.

There are a number of sources of relatively intense electromagnetic fields in the everyday environment. Some of these sources are obvious and easy to avoid. Others are not apparent and exposure is unavoidable. However, we believe that by following the warnings listed below, your risk to EMI will be minimized.

The sources of radiated EMI can be broadly classified into three types

1. Hand-held portable transceivers (transmitters-receivers) with the antenna mounted directly on the transmitting unit. Examples include: citizens band (CB) radios, "walkie talkie," security, fire, and police transceivers, mobile phones, and other personal communication devices.



- Some mobile phones and similar devices transmit signals while they are ON, even when not being used
-

2. Medium-range mobile transceivers, such as those used in police cars, fire trucks, ambulances, and taxis. These usually have the antenna mounted on the outside of the vehicle; and
3. Long-range transmitters and transceivers such as commercial broadcast transmitters (radio and TV broadcast antenna towers) and amateur (HAM) radios.

2. IMPORTANT PRECAUTIONS



- Other types of hand-held devices, such as cordless phones, laptop computers, AM/FM radios, TV sets, CD players, and cassette players, and small appliances, such as electric shavers and hair dryers, so far as we know, are not likely to cause EMI problems to your mobility scooter.
-

MOBILITY SCOOTER ELECTROMAGNETIC INTERFERENCE:

Because EM energy rapidly becomes more intense as one move closer to the transmitting antenna (source), the EM fields from hand-held radio wave sources (transceivers) are of special concern. It is possible to unintentionally bring high levels of EM energy very close to the motorized scooter control system while using these devices. This can affect mobility scooter movement and braking. Therefore, the warnings listed below are recommended to prevent possible interference with the control system of the mobility scooter.

Warnings:

Electromagnetic interference (EMI) from sources such as radio and TV stations, amateur radio (HAM) transmitters, two-way radios, and mobile phones can affect mobility scooters. Following the warnings listed below should reduce the chance of unintended brake release or mobility scooter movement which could result in serious injury.

1. Do not operate hand-held transceivers (transmitters-receivers), such as citizens band (CB) radios, or turn ON personal communication devices, such as mobile phones, while the mobility scooter is turned ON;
 2. Be aware of nearby transmitters, such as radio or TV stations, and try to avoid coming close to them;
 3. If unintended movement or brake release occurs, turn the mobility scooter OFF as soon as it is safe;
 4. Be aware that adding accessories or components, or modifying the mobility scooter, may make it more susceptible to EMI; and
-



- There is no easy way to evaluate their effect on the overall immunity of the mobility scooter.
-

5. Report all incidents of unintended movement or brake release to the distributor listed on the inside front cover of this manual. Note whether there is a source of EMI nearby.




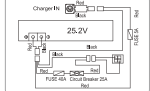













Important Information:

1. 20 volts per meter (V/m) is a generally achievable and useful immunity level against EMI (as of May 1994). The higher the level, the greater the protection.
2. The immunity level of this product is at least 20/Vm.

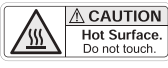2. IMPORTANT PRECAUTIONS

2-7 SAFETY WARNING AND INSTRUCTION LABELS





<p>1</p> 	<p>Warning Sticker</p> <ol style="list-style-type: none"> 1. Please read the instruction booklet carefully before using your scooter. 2. Keep the scooter properly maintained. 3. Do not drive the scooter on slippery surfaces or on slopes over 6°. 4. Do not drive on highway, crowded roads, or unfamiliar areas. 5. Do not turn at high speed in either forward or reverse. 6. Do not wash with water or leave scooter in humid environment since water can damage the electronic parts. 7. Always re-engage the emergency freewheel drive after use 8. Weight restriction of the device : 255 lbs 9. Do not operate on staircases or escalators. 		
<p>2</p> 	<p>Scooter information label</p> <p>States the date of manufacture and serial number of the scooter. This product has been supplied from an environmentally aware manufacturer. It may contain substances that could be harmful to the environment. Recycling must be carried out in accordance with the respective national legal provisions.</p>		
<p>3</p> 	<p>Power warning label</p>	<p>11</p> 	<p>Wiring diagram for batteries</p>
<p>4</p> 	<p>Instruction of battery lock nut (circuit breaker) label</p>	<p>12</p> 	<p>Seat position warning plate</p> <p>Do not get on the scooter if the seat is not fixed. Push the seat forward until the seat locks to its position (green), you will hear a click sound.</p>
<p>5</p> 	<p>Battery charging warning label</p>	<p>13</p>  <p>Automatic Version Only</p>	<p>Folding / unfolding button instruction label :</p> <p>Press folding or unfolding button (either by pressing button on scooter or by remote) once for 3 seconds to wake up or turn off the folding / unfolding function.</p>
<p>6</p> 	<p>European representative label</p>	<p>14</p> 	<p>Safety instruction label of Lithium battery :</p> <p>Read carefully to prevent improper operation.</p>
<p>7</p> 	<p>N-D Lever Adjustment label</p> <p>Instructs freewheel mode operation. (see section 4-4)</p>	<p>15</p> 	<p>Fixation point reminder</p> <p>Sit only when the seat slide is fixed to green point</p>
<p>8</p> 	<p>Heat caution warning label</p> <p>Be aware that the motor may normally become hot during operating. Please avoid skin contact.</p>	<p>16</p> 	<p>Seat position indication label</p> <p>Seat latch is locked when the seat positioned in green. Seat latch is released when the seat positioned in red. Do not sit if seat is positioned in red</p>
<p>9</p> 	<p>Pinch point warning label</p> <p>Keep hands away</p>	<p>17</p> 	<p>Pinch point warning label</p> <p>Keep hands away</p>
<p>10</p> 	<p>Charge port cap sticker</p>		

3. IDENTIFICATION OF PARTS



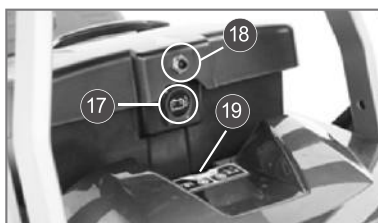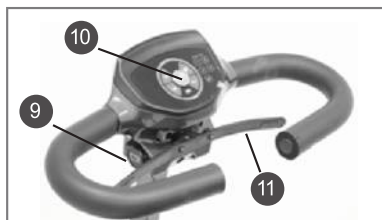Figure 4



No.	ITEM	No.	ITEM
01	Steering	13	Seat Backrest
02	Tiller Height Adjustment	14	Seat Knob
03	Head Light	15	Battery Box
04	Front Wheels	16	Battery Box Handle
05	Footplate Cover	17	Charging Port (battery box)
06	Seat Springs	18	Circuit Breaker
07	Supporting Wheels	19	Folding / unfolding light indicator / buttons
08	Rear Wheels	20	Tail light / Brake light
09	Charging Port (upper)	21	N-D Lever
10	Control Panel / Power Button	22	Anti-Tippers
11	Throttle Lever	23	Rear Bumper
12	Seat		

*Above is only for reference, actual product may vary.

4. INTRODUCTION & OPERATION OF SCOOTER

4.1 IDENTIFICATION & OPERATION OF CONTROL PANEL

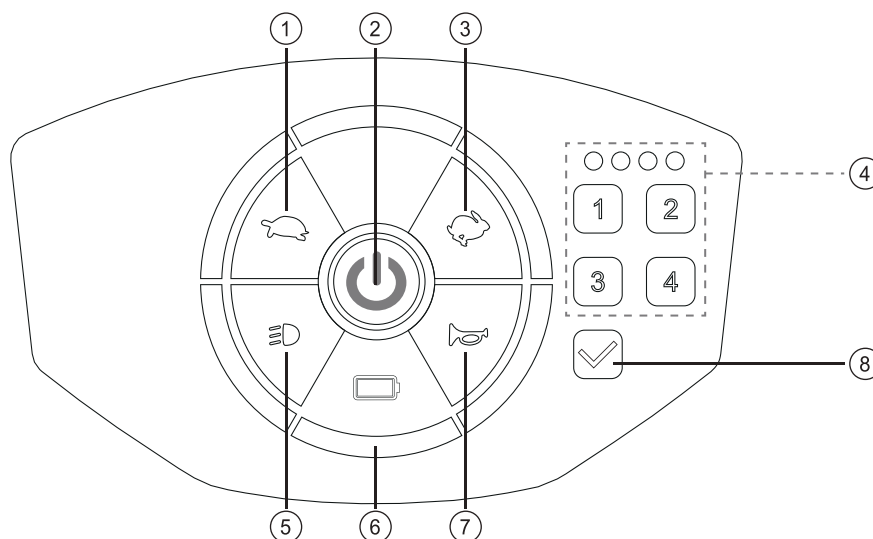


Figure 5 - Control Panel

No	Item	Function
1	Speed Deceleration	Press to decrease speed
2	Power Switch	Press to turn on / off the power
3	Speed Acceleration	Press to increase speed
4	Digital Lock (number pad)	Set a password to lock the scooter
5	Head-light / Tail-light Switch	Press to turn on headlight and tail-light
6	Battery Gauge / Error Light Indicator	Showing battery status or errors
7	Horn	Press to sound warning tone when necessary
8	Digital Lock (enter key)	Press Enter when password is confirmed



- Please keep the panel clean, away from water, oil, and dust as the RF reading range may be reduced or malfunction may result.
- When operating the panel, always press on the center of the button to ensure the button press is detected.
- Do not place any objects on the panel or it may disrupt the touch screen or cause damage to the scooter.
- Battery / Error Indicator will be flashing if there is error, please stop operating the scooter if indicator is flashing (see section 8).
- If the power button continues lightening after power off, please release the throttle control lever and remove the battery pack. Please contact your dealer.

4. INTRODUCTION & OPERATION OF SCOOTER

POWER SWITCH:

- Press the button to turn on/off the power
- Power standby (locked) - Red light
- Power on (unlocked) - Green light
- Power off - Light off

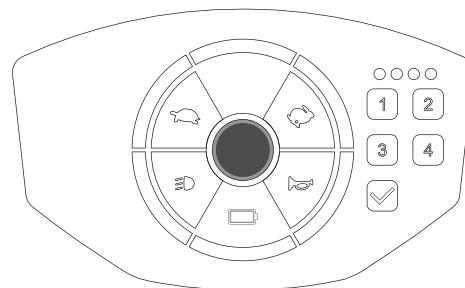


Figure 6



Auto-power off :

Power will automatically be off in 30 seconds if no activity is detected during power standby (locked, red light) period to prevent tap error.

DIGITAL UNLOCK YOUR SCOOTER:

- Press the power switch, the power (red light) will be on standby for you to unlock.
- To unlock, enter your password on the number pad (D), light indicator (C) will light up each time when the number is entered. Press (E) on the pad when password entering is finished. (Figure 7)
- The power light will turn green when the scooter is successfully unlocked..

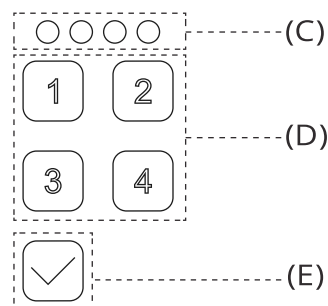


Figure 7



Factory Default Password : 1234

If there is enter error, press the (E) key or power switch button again to re-enter the password.

PASSWORD SETTING:

- Once scooter is unlocked (power green light on), press and hold the (E) key for 3 seconds, the pad light will be flashing three times.
- Enter your password on number pad (D), the password can be 1-4 digits number, digits can be repeated. Light indicator (C) will light up each time when the number is entered.
- Press (E) on the pad when password entering is finished, pad light will be flashing three times.



If there is enter error, press power switch button again or wait for 30 seconds to re-enter the password. Scooter must be unlocked before reset the password. If you forget the password, please unlock the scooter by RF sensor key.

4. INTRODUCTION & OPERATION OF SCOOTER

RF SENSOR UNLOCK YOUR SCOOTER:

- Press the power switch, the power (red light) will be on and standby for you to unlock.
- To unlock, place the RF sensor tag to the panel (spot near power switch for best reading result) (Figure 8)
- The power light will turn green when the scooter is successfully unlocked.

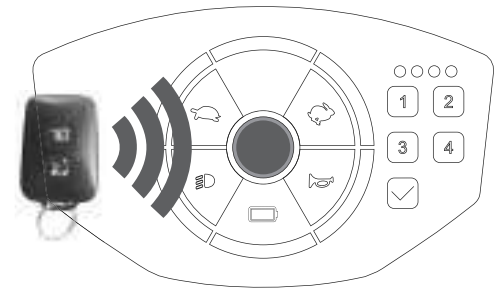


Figure 8



- The sensor transmit is up to 20mm, however, please keep the distance between RF sensor and panel as close as possible.
- Do not disassemble the RF remote, the inner sensor chip may be damaged.
- Do not turn off the power while unlocking the scooter, no matter it's Digital Unlock or RF Sensor Unlock. User will need to turn the power back on and re-unlock again if unlocking is interrupted.
- Scooter cannot be operated unless it's unlocked (green light on).

PANEL BACK-LIGHT ADJUSTMENT (FIVE LEVELS):

- Press the power switch to get the scooter on power standby (red light on), then press and hold both acceleration [🏍️] / deceleration [🛑] buttons for 2 seconds. The panel light will be all on when scooter is ready for back-light setting. (Figure 9)
- Press [🏍️] to increase level of back-light. Press [🛑] to decrease level of back-light.
- Press the power switch to save the back-light setting
- Panel back-light can be adjusted under power standby mode (red light), user doesn't have to unlock the scooter.

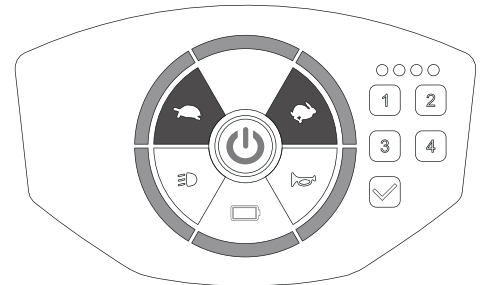


Figure 9

BATTERY / ERROR LIGHT INDICATOR:

- The light will be showing the battery status or errors of scooter. If scooter has errors, the light indicator will be flashing (see section 8). This light indicator is also the indication of battery status (see section 6-3). (Figure 10)

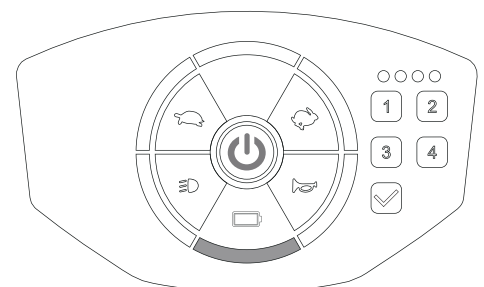


Figure 10

4. INTRODUCTION & OPERATION OF SCOOTER

SPEED ADJUSTMENT:

- Press acceleration [🐢] / deceleration [🐌] button to increase / decrease speed.
- There are 3 different speeds of adjustment fastest 3 > 2 > 1 slowest. (Figure 11)

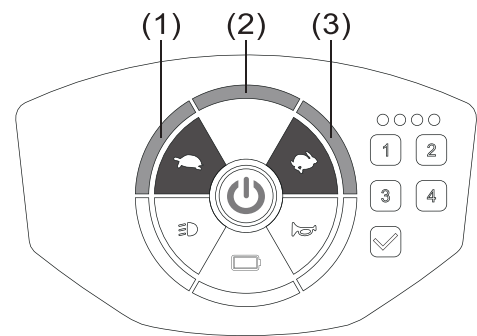


Figure 11



- Do not set the speed at maximum when operating the scooter.
- Adjust the speed at maximum when inclining a slope, and adjust to minimum speed when declining.

HEAD-LIGHT / TAIL-LIGHT:

- Press button [💡] once to switch on the head-light and tail - light, the [💡] will be lighted on the panel.
- Switch off by pressing the button again, the [💡] light will be off on the panel (figure 12)

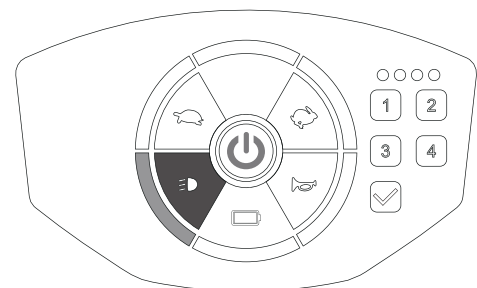


Figure 12



- Tail-light will be on anyway when the brake is engaged and when the scooter speed is reduced.

HORN:

- Press horn button to sound the buzzer. Buzzer sound and light will last until you release the horn button (Figure 13)

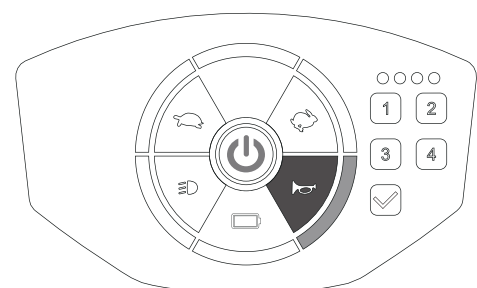


Figure 13

4.2 OPERATION OF MAIN SWITCH:

MAIN SWITCH (TOUCH SCREEN VERSION):

- Press power button to turn on/off (Figure 14) (see section 4-1 for details)



Figure 14

4. INTRODUCTION & OPERATION OF SCOOTER

MAIN SWITCH (KEY VERSION):

- Turn the key clockwise to turn on, counterclockwise to turn off (Figure 15)

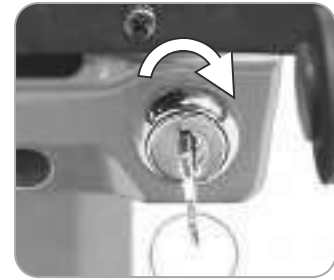


Figure 15

SLEEP MODE:

- Scooter will go into sleep mode with one long beep warning sound if no throttle activity is detected for programmable time period. (Default programmable time is 10 mins.)
- When scooter is in sleep mode, all lights on the control panel will be off, and the scooter will not respond to commands.
- To wake up the scooter, turn the power (key) switch off and then on again.

4.3 OPERATION OF THROTTLE:

THROTTLE:

- Pull the right throttle (F) to move scooter forward. Pull the left throttle (G) to move scooter backward. This can be reversed if required by local dealer. (Figure 16)
- These are also your accelerator. The further you depress them, the faster you go. (Subject to the level of speed adjustment)
- Releasing both left/right throttles will automatically operate the brakes to slow down and stop.

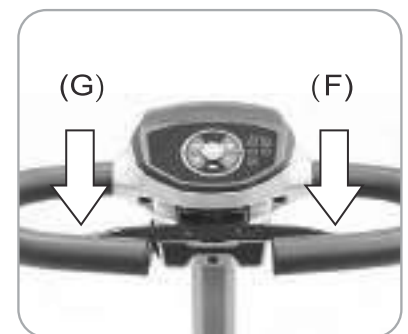


Figure 16



- Do not pull both right and left throttle at the same time; you might not be able to control the scooter.
- Turning the scooter OFF whilst driving will bring the scooter to an abrupt stop and danger.

4. INTRODUCTION & OPERATION OF SCOOTER

4.4 OPERATION OF FREE-WHEELING (N-D LEVER)

FREE-WHEELING (N-D LEVER):

- When lever is in Drive (D) position, the scooter can be driven. Normal position is in D (Figure 17)
- When lever is Neutral (N) position, the scooter can be moved manually without power (Figure 18)



Figure 17



Figure 18



- Freewheel operation is only recommended on flat surfaces. Never leave your scooter on gradient with its motors disengaged.
- When adjusting N-D lever, do not sit on the scooter.
- It's not able to drive the scooter when the lever is in Neutral. You must turn off the scooter first and switch to D lever, then turn on and drive scooter.

4.5 TILLER HEIGHT ADJUSTMENT

- Loosen (counterclockwise) the screw knob to adjust the tiller height, tighten (clockwise) screw knob again to lock at a desired position.
- There is a stop point on tiller, the height cannot exceed the stop point (approx. 20cm) (Figure 19)

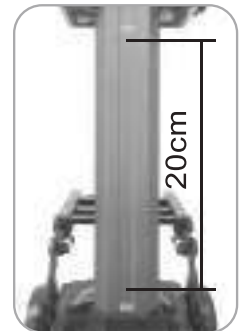


Figure 19

4.6 SEAT ADJUSTMENT / DISMANTLE

- Please refer to section 5 for scooter folding / unfolding.
- When scooter is being folding, the seat latch will automatically release and seat will move backwards.
- After unfolding, the seat shall automatically move forwards. However, please always double check and push the seat to its locked position) manually if not.
- Do not get on the scooter if the seat is positioned in red point(H).Get on the scooter only if the seat is locked in its position green point (Fig.20)



Figure 20

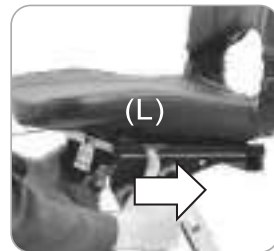
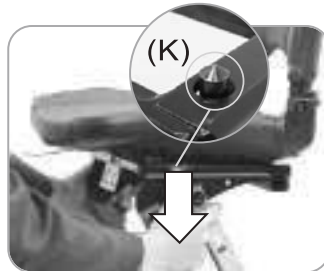
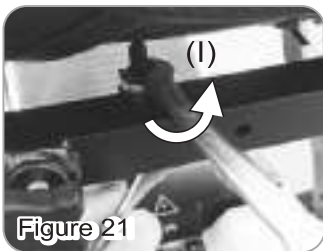
4. INTRODUCTION & OPERATION OF SCOOTER

Seat Dismantle (Figure 21.22)

1. Loosen the screws (I) by wrench and hex key.
2. Press the springs (J) to remove from the seat slide bracket.
3. Pull down and hold the seat springs to release the seat latch (K).
4. Push the seat bracket (L) by thumbs towards the back of scooters to release the seat from latch.
5. When the seat is released from latch, grab the seat (L1) and push it towards the back of scooter to remove the seat.

※When putting the seat back, insert the seat back to seat slide bracket. Push the seat towards to the front of scooter (green point), fasten the screws (Fig.21) and insert the springs back to seat slide bracket (J).

※Please make sure to follow the direction of Step 2 (Fig. 22). shown when putting springs (J) back to the bracket. Otherwise, it may cause damage.



- Do not conduct any adjustment (tiller / seat / armrest) while driving the scooter.
- Do not get on the scooter until the seat firmly locks in the position (green point).
- Except for removing seat, do not operate the seat sides spring if no necessary.
- When pulling the seat springs, please grab and hold down both sides springs, this helps user to easily release the seat latches. The latches may not be release if the holding point of seat springs is incorrect.

5. FOLDING / UNFOLDING YOUR SCOOTER

Folding / unfolding your scooter only when the scooter is under below circumstances:

- When the scooter is placed on a flat ground and in lay down position.
- When the scooter is power off and on standby for folding / unfolding, the folding / unfolding light indicator (M) on scooter will be orange. (Figure 23)

Folding / Unfolding Buttons:

- Buttons on scooter (Figure 23)
 - Press folding button (N) once to auto-fold; please make sure the scooter is fully expanded before folding.
 - Press unfolding button (O) once to auto-unfold, please make sure the scooter is fully folded before unfolding.
 - Press (N) / (O) again to fold / unfold if folding / unfolding process is interrupted.
- Buttons on remote (Figure 24)
 - Press folding button (P) once to auto-fold, please make sure the scooter is fully expanded before folding.
 - Press unfolding button (Q) to auto-unfold, please make sure the scooter is fully folded before unfolding.
 - Press (P) / (Q) again to fold / unfold if folding / unfolding process is interrupted.

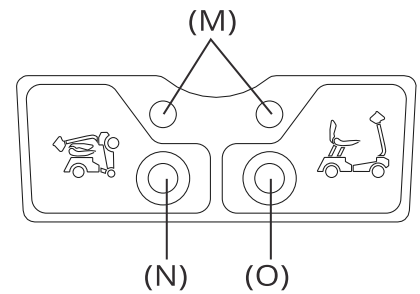


Figure 23

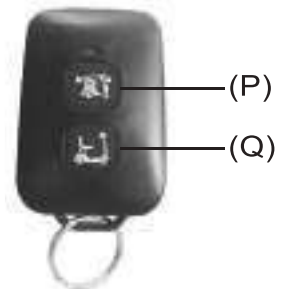


Figure 24



- The remote is not designed to be water-proof, please keep it away from humid environment.
- Best transmit range of remote is within 10m. (The reading range will be affected due to the environment, barrier, and weather).
- The scooter may not be folded / unfolded if the battery is low, please re-charge the scooter by upper tiller charging port immediately and operate the folding / unfolding later when the scooter is fully charged. (see section 6-4)
- Keep the remote carefully, it's highly recommended to chain with RF sensor key.
- Please contact dealer if the key is lost.

5. FOLDING / UNFOLDING YOUR SCOOTER

How to wake up the scooter:

- Scooter will go into sleep if no throttle activity is detected for programmable time period. (Default programmable time is 3 mins.)
- Please wake up the scooter before folding / unfolding. There are 3 ways to get the scooter back on standby for folding/unfolding:
 1. Press power switch button and turn off the power.
 2. Press and hold the folding button (N) or unfolding button (O) on scooter for 3 seconds to wake up the folding / unfolding function. The light indicator (M) will be orange. To turn off the folding / unfolding function, press and hold the folding button (N) or unfolding button (O) on scooter for 3 seconds again.
 3. Press and hold the folding button (P) or unfolding button (Q) on remote for 5 seconds to wake up the folding / unfolding function. To turn off the folding / unfolding function, press and hold the folding button (N) or unfolding button (O) on remote for 5 seconds again.

How to fold the scooter (Figure 25):

1. Turn off the power
2. Flip down the seatback
3. Press folding button on the remote or on scooter once to auto-fold the scooter. The seat will automatically moving backwards during folding.
4. The seat latches will automatically release and seat will move backwards.



Figure 25

5. FOLDING / UNFOLDING YOUR SCOOTER

How to unfold the scooter (Figure 26):

1. Place the scooter to a flat ground and lay down the scooter
2. Make sure the scooter is power off, however is awake for unfolding (the light indicator M will be orange).
3. Press unfolding button on the remote or scooter once to auto-unfold the scooter.
4. When folding and unfolding, the seat shall automatically shift to its position, however, please double check and push seat to its locked position (green) manually if not.
5. Flip up the seatback and double check if lock is all in its position before you get on the scooter.



Figure 26



- Please turn off the power before folding / unfolding, the scooter can only be folding / unfolding when it is power off, the folding / unfolding light indicator will be orange.
- There is protection for tapping error; scooter cannot be unfolding when it is in standing position.
- Actions will stop if any foreign objects are detected during folding / unfolding process. Please remove the objects and conduct the actions again.
- When operating the remote, please stand beside the scooter for best receiving signal.
- When operating the folding / unfolding buttons on scooter, please reach the buttons from the back of scooter. User may get pinched if reaching buttons from other sides of scooter.
- The scooter may not be folded / unfolded if the battery is low, please re-charge the scooter by upper tiller charging port immediately and operate the folding / unfolding later when the scooter is fully charged. (see section 6-4)
- Check and replace the battery for remote control if the scooter is not able to fold /unfold.
- After completion of unfolding the scooter, it is normal that the scooter foot area will slightly sink while seating for the first time.

6. BATTERY GAUGE & CHARGING THE BATTERIES

6.1 DISASSEMBLING BATTERY BOX

- Each battery box contains one lithium battery or (2) lead-acid batteries (optional).
- The battery box is designed for a better stability. Please follow the instruction to disassemble the battery box (Figure 27) . Push the battery backwards (1.) and pull it upwards (2.).
- It's highly recommended not to remove battery box out if no necessary. Always keep battery box installed on scooter.
- For any reason to remove the battery box, please make sure the connectors are in sequence and in position when put it back to scooter.

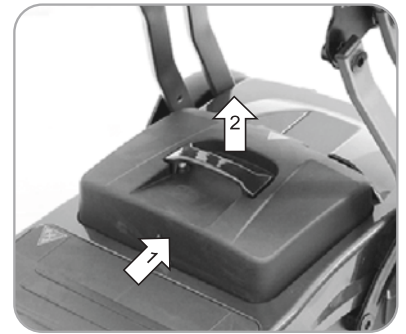


Figure 27



- When putting battery box back to scooter, gently press battery box downward to make sure it is firmly connected to the scooter.
- Do not disassemble battery box and open sealed parts by yourself to prevent electric shock and damage the battery.
- Maintenance & repair must only be carried out by a competent engineer or authorized dealer. Please contact dealer if you need to replace internal battery.

6.2 CAUTIONS FOR BATTERY & CHARGER

Lithium Battery:

- Only certain type and voltage of battery can be used on the scooter. Your scooter may supply with standard lithium battery or optional lead-acid battery.
- Scooter can be charged when it's folded or unfolded
- Only use the lithium battery charger (under 4 Amps) supplied with your scooter for lithium battery.
- Do not switch off, unplug or interrupt the recharge cycle until the charging cycle has completed (charger indicator turns green)
- If storing your scooter for a prolonged period (1 month or more), storage the batteries at 30%-50% state of charge at least. Batteries must be recharged once every three months at least if storing for more than six months. On returning, make sure to charge them again before using scooter.
- Store the scooter in a dry, well-ventilated environment. Do not store your scooter where it will be exposed to source of direct heat, damp, oil, acid, alkaline, or where Ozone could be possibly generated. All of the above will minimize scooter's lifetime.

6. BATTERY GAUGE & CHARGING THE BATTERIES

- Recommended range of storage : Temperature : 20~30°C, Humidity : 25% ~ 85%, out of this range might affect the performance of the charger and battery.
- Rough handling operations, throwing, rolling, crushing, shaking, overly strong impact and extrusion to batteries should be avoided, to prevent battery casing broken or internal structure damage.
- It's highly recommended to have your scooter inspected for battery every six month despite how often you use the scooter. Please contact your dealer.
- Batteries are subject to disposal and recycling regulations that vary by country and region. Always check and follow your applicable regulations before disposing of any battery.

24V/12A Lead-Acid Battery (optional):

- Only certain type and voltage of battery can be used on the scooter. Your scooter may supply with standard lithium battery or optional lead-acid battery.
- Only use the lead-acid battery charger (under 8 amps) supplied with your scooter for lead-acid battery
- Adding acid electrolyte to battery is prohibited
- Batteries will only give maximum performance after scooter has been used and batteries have been recharged several times.
- The minimum time needed to recharge will vary depends on the depletion of the batteries. Excessive or short period charging will be detrimental to battery life.
- Charge the batteries after each trip. Do not turn the key on while scooter is in charging.
- Please charge the scooter in 0°C ~ 40°C, out of this range might affect the performance of the charger and battery.
- If storing your scooter for a prolonged period, batteries must be recharged once every three months at least
- For longest life, your batteries should be recharged regularly. Please recharge the batteries before they run down to 20%.
- It's highly recommended to have your scooter inspected for battery every six month despite how often you use the scooter. Please contact your dealer.

Chargers(lithium / lead-acid):

- Do not leave charger plugged into your scooter with charger switched off as this may discharge your batteries. Always unplug the cables when finishing charging.
- Inspect the battery charger before each use, make sure connectors are dry and clean. Do not attempt to use an extension cord to plug in your battery charger.
- Be aware that the battery charger case may become hot during charging. Please avoid skin contact.

6. BATTERY GAUGE & CHARGING THE BATTERIES

- Do not cover charger by any objects. Use the charger in a well-ventilated area. Do not smoke as explosive gases may be generated while charging the batteries
- Keep the charger stored in a well-ventilated area, do not drop the charger.



- Lead-acid battery charger cannot be alternatively used on lithium battery. Please only use specific lithium battery charger for your lithium battery. Lead-acid battery charger can be used only if your scooter comes with lead-acid batteries (optional). (Please see section 9-1 for details)
 - If you want to change to optional lead-acid batteries, please contact dealer to reset the battery system such as charging port and controller setting.
 - Please keep the above cautions in mind. Any battery / charger faults due to the unauthorized maintenance, dismantle, misuse or accidental damage is not covered by the manufacturer's warranty.
-

6.3 BATTERY GAUGE

There is an indicator shows battery capacity status (Figure 28):

- **Touch screen version : 6 segments**
- When battery gauge light trends to right (green), it indicates sufficient power capacity.
- When battery gauge light trends to left (red), it indicates insufficient power capacity.



Figure 28



- Check battery capacity before driving. Recharge the battery immediately when the red light is flashing as the remaining power can only keep scooter going for buffering 2 km or less
 - Displayed segments on panel are for reference only. Please refer to LED indicator of charger for more accurate charging status.
 - The scooter may not be able to fold/unfold when the battery capacity status is less than 3 segments, please recharge the battery immediately.
-

6. BATTERY GAUGE & CHARGING THE BATTERIES

6.4 OPERATING THE CHARGER:

Batteries must be charged before using the scooter for the first time and should be recharged after each day use.

Your scooter may come with standard lithium 24V/2A charger, or optional bigger 24V/4A charger if you upgrade the lithium battery to 17.4Ah . Your scooter may come with 24V/2A lead-acid charger if you optionally choose lead-acid battery for your scooter. Do not alternatively use lead-acid charger to charge lithium battery, they are not compatible and may cause damages and dangers if doing so.

Lithium Battery Charger:

- Output Voltage: 100V~240V (Figure 29)
 1. Charger output cable, 4pins ☹ connector (R)
 2. AC power cord (S)
 3. LED Indicator (T)
 - RED LIGHT ON : Pre-charging
 - ORANGE LIGHT ON : Normal Charging
 - GREEN LIGHT ON: Fully Charged.
 4. 4A takes 5 hours for a full charge. Remove the AC plug first after charging.

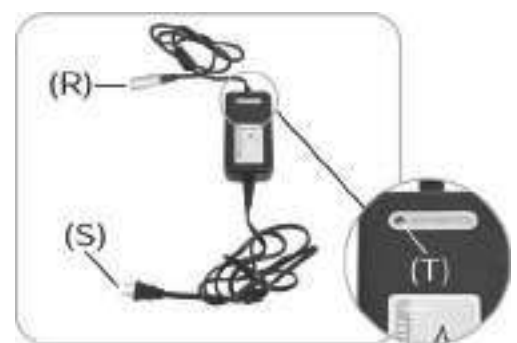


Figure 29

Lead-acid Battery Charger:

- Output Voltage: 100V~240V (Figure 30)
 1. Charger output cable, 3pins ☹ connector (U)
 2. AC power cord (V)
 3. LED Indicator(W)
 - RED LIGHT ON : Pre-charging
 - ORANGE LIGHT ON : Normal Charging
 - GREEN LIGHT ON: Fully Charged.
 4. takes 7 hours for a full charge. Remove the AC plug first after charging.

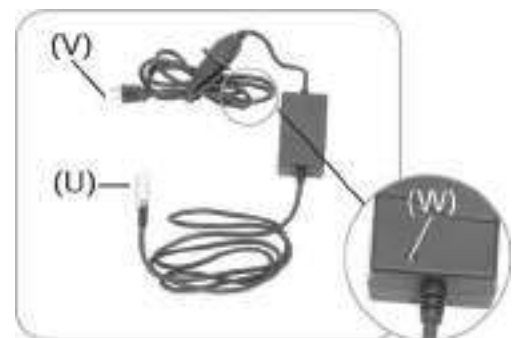


Figure 30



- Each country may supply different charger. The charging procedure may be different from below. Chargers for lithium battery and lead-acid battery (optional) are different.
- The charger LED light should indicate red when pre-charging. If not, please contact your authorized dealer.
- Inspect the battery charger before each use, make sure connectors are dry and clean. Do not attempt to use an extension cord to plug in your battery charger.
- Do not use other brand of charger to charge the scooter, only use the charger equipped with originally

6. BATTERY GAUGE & CHARGING THE BATTERIES



- The battery charger is equipped with cooling fan. If the fan does not appear to be working correctly, please turn off the charger immediately as it may be overheated. Please contact your authorized scooter dealer.
- Do not leave charger plugged into your scooter with charger switched off as this may discharge your batteries. Always unplug the cables when finishing charging.
- Be aware that the battery charger case may become hot during charging. Please avoid skin contact.
- Keep the charger stored in a well-ventilated area, do not drop the charger.

Charging Operation:

There are two charging ports, one is on upper of scooter, and the other is on battery box. Upper Charge Port (Figure 31)

- Open the charge port cap
 1. Insert the charger output cable into the charge port
 2. Plug the AC cord into a standard electrical outlet on the wall
 3. Check the charger LED indicator; make sure it's in normal charging. During charging, LED will indicate orange light.
 4. When unplug, please make sure to remove the AC cord from the wall first, and then charger output cable from the charge port.

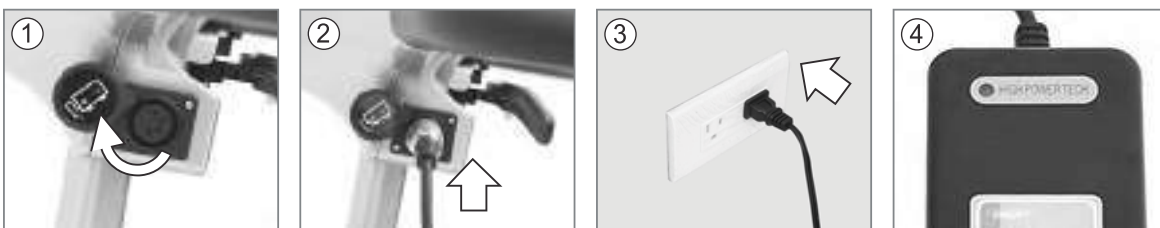


Figure 31

6. BATTERY GAUGE & CHARGING THE BATTERIES

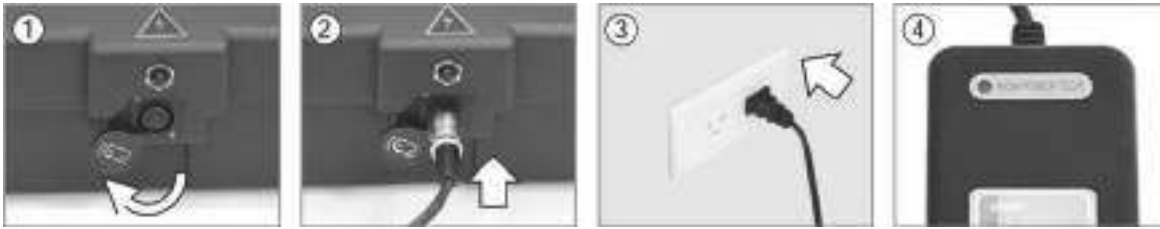


Figure 32

- Battery Box Charge Port (Figure 32)

1. Remove the battery box from the scooter and place it on a flat ground.
2. Open the charge port cap and connect the charger output cable into the charge port
3. Plug the AC cord of charger into a standard electrical outlet
4. Check the charger LED indicator; make sure it's in normal charging. During charging, LED will indicate orange light.
5. When unplug, please make sure to remove the AC cord from the wall first, and then charger output cable from the charge port.
6. Put the cap back to cover charging port before install the battery box back to scooter.



- Be sure the scooter key (power) is in the OFF position when charging.
- The time needed to recharge will vary depending on the depletion of the batteries. Charging for longer than necessary will not harm the batteries. It cannot be over-discharged.
- When remove / install the battery box, please make sure the connectors are in sequence and in position to avoid electric damage.
- Please keep the above cautions in mind. Any battery / charger faults due to the unauthorized maintenance, dismantle, misuse or accidental damage is not covered by the manufacturer's warranty
- Do not alternatively use lead-acid charger to charge lithium battery, they are not compatible and may cause damages and dangers if doing so.

7. CARE AND MAINTENANCE

7.1 INSPECTION CHECKS

Daily Check :

Please always check your scooter before you start your every journey.

Check point	Inspection	What to do if the inspection is failed
N-D leaver	Check for correct function.	Contact your dealer.
Horn	Check for correct function.	Contact your dealer.
Throttles	Pull the wigwag to test the scooter movement. Pull the wigwag a little bit and release it to test if brake works.	Contact your dealer, the electromagnetic brake should be clean by brush regularly.
Remote / folding-unfolding buttons on scooter	Check if Remote control for folding / unfolding works. Check if folding / unfolding buttons on scooter works.	<ul style="list-style-type: none">• Replace the battery for remote control.• Contact your dealer.
RF remote control	Check for correct function.	Change battery (3Vx1pcs) or contact your dealer.
Folding / unfolding Release lever (manual version)	Check if folding / unfolding function is normal.	Contact your dealer.
Battery Gauge	Check if the battery gauge is displayed and whether it is at low power.	<ul style="list-style-type: none">• Contact your dealer if battery gauge is not working.• Recharge the battery immediately if low.
Battery box	Check if there is cracks or damage, and if it's firmly installed to the scooter.	Contact your dealer.
Lighting	Check if all lights, such as head lights, rear lights, and turn signal are functioning correctly.	Contact your dealer.
Armrests	Check if the parts are clean and firmly tighten to the scooter and do not wobble. Tighten the screw knob that holds the armrest.	Contact your dealer.
Wheels / Tyres	Check if the tyre has cracks, deformation, and check that wheels rotate without wobbling. Tyre tread depth is not less than 0.5mm. No foreign objects in tyres.	Contact your dealer.

7. CARE AND MAINTENANCE

Check point	Inspection	What to do if the inspection is failed
Motor	Check for correct function, no noisy sound.	Contact your dealer.
Battery Charger	Check if the charger is functioning correctly and the batteries are charged.	Contact your dealer.

MONTHLY CHECK

Check point	Inspection	What to do if the inspection is failed
Seat / Upholstery	Check for movement and seat lock and if upholstery is worn.	Fasten the screws if it's loosen Contact your dealer if upholstery is worn.
Tiller fixation	Check for all fasteners.	Contact your dealer.

EVERY SIX MONTHS CHECK

Check point	Inspection	What to do if the inspection is failed
Electronics	Check if all the battery cables and connectors are firmly tighten to the scooter.	Contact your dealer.

Please take your scooter to a qualified technician for a thorough inspection and servicing every six months.

7.2 CLEANING AND MAINTENANCE

- User should inspect the scooter regularly to keep scooter in good running order.
- User should inspect the scooter more often especially when the scooter has possibly been riding on the bad road condition such as, grass, sand, bumps, and potholes.
- Please turn off the scooter and remove the key and make sure the lever is in D position before inspection to avoid roll-away.
- Please avoid any tapping error on remote or scooter for folding/unfolding when carrying the maintenance work.
- Do not use any abrasive or scouring liquids for cleaning. Only use a damp cloth and gentle detergent.
- Do not use hose pipe or splash water directly onto the scooter as this may cause damage to electronics.

7. CARE AND MAINTENANCE

Seat Upholster:

Only use damp cloth and a mild soap to wipe the seat. Do not use abrasive cleaners as this will damage the seat. Do not sit until the seat dries out.

Frame / Shroud:

Only use damp cloth to wipe. Do not use abrasive cleaners as this will damage the parts. When scooter dries out, the shroud can be applied with a light of wax to keep its gloss.

Tyre:

Only use damp cloth to wipe. Do not use abrasive cleaners as this will damage the parts. User should inspect the tyres frequently for damage, the presence of foreign bodies, unusual wear and sufficient tread depth. If replacement of tyres are needed, please contact the nearest dealer.

The following areas require periodic inspection: Tread depth drops below 0.5 mm.

Storage:

- If store the scooter in long time (one week or more), please disconnect the battery terminals and do not connect battery terminal [+] [-] with any metals to avoid danger.
- When store scooter, make sure the power is off, seat lock is fixed, and N-D lever is in D position to engage brake to avoid roll-away.
- Please store the scooter in a dry location. Do not store your scooter where it will be exposed to source of direct heat, damp, oil, acid, alkaline, or where Ozone could be possibly generated. All of the above will minimize scooter / tyre cycle and shorten its lifetime. Recommended range of storage for the scooter is as below :
- Temperature (with lithium battery) : Less than 1 year : -20°C ~ 20°C
 - Less than 3 months : -20°C ~ 40°C
 - Less than 1 month : -20°C ~ 50°C
- Temperature (with lead-acid battery) : -30°C ~ 50°C
- The highest temperature during battery transportation should be lower than 60°C
- If the storage environment is out of the range, we suggest you to take out the battery box from scooter and store it respectively.

All maintenance and repair of scooter MUST be done by an authorized dealer. ANY unauthorized adjustments, repair or service to the scooter may result injury to yourself or damage to scooter.

Here are some suggestions about solving problems you may have with your scooter. There is a self-diagnostic warning light on the control panel. To check the self-diagnostic warning light, turn on the key (power) switch (see section 4-2) and count the number of blinks on the warning light.

8. TROUBLE SHOOTING

Touch Screen version:

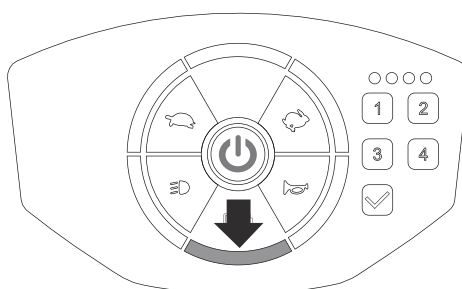


Figure 33 Normal : Solid light
Abnormal : flashing light

Self-Diagnostic Warning Light

Flash	Description	Initial check points
1	Battery Low	The batteries are running low. • Recharge the batteries.
2	Low Battery Fault	The batteries have run out of charge. • Recharge the batteries. • Check the battery and associated connections and wiring.
3	High Battery Fault	Battery voltage is too high. This may occur if overcharged and/or traveling down a long slope. • If traveling down a slope, reduce your speed to minimize the amount of regenerative charging.
4	Current limit time-out or controller overheat	The motor has been exceeding its maximum current rating for too long. • The scooter may have stalled. Turn the controller off, leave for a few minutes and turn back on again. • The motor may be faulty. Check the motor and associated connections and wiring.
5	Park Brake	Either a park brake release switch is active or the park brake is faulty. • Check the park brake and associated connections and wiring. • Ensure any associated switches are in their correct positions. • Switch to D (drive) position. Turn off the power and turn on again.
6	Drive Inhibit	Either a stop function is active or charger inhibits or OONAPU condition has occurred. • Release the stop condition (seat raised etc.) • Disconnect the battery charger. • Ensure the throttle is in Neutral when turning the controller on. • The throttle may require re-calibration. • The folding/unfolding procedure is uncompleted.
7	Speed Pot	The throttle, speed limit pot. SRW or their associated wiring may be faulty. • Check the throttle and speed pot and associated connections and wiring.
8	Motor Voltage	The motor or its associated wiring is faulty. • Check the motor and associated connections and wiring.
9	Other Error	The controller may have an internal fault. • Check all connections and wiring.

- If self-diagnostic warning lights with 1 / 2 / 4 / 5 times, please refer to above for solutions (turn off and restart again or charge the battery). If none of above correct problem, contact your authorized dealer.
- If self-diagnostic warning lights with 3 / 6 / 7 / 8 / 9 times, please contact your authorized dealer directly.

8. TROUBLE SHOOTING

Folding / Unfolding Errors:

Here are some suggestions about solving problems you may have when folding / unfolding your scooter. The scooter will stop the folding / unfolding actions if errors are detected. Please refer to below table and correct the problems.



Figure 35



Figure 36

Flash	Description	What to do
1	Power On	Please turn off the power of scooter (see section 4-2) and restart the folding / unfolding procedure.
2	Battery Low	The batteries are running low. Recharge the batteries and restart the folding/unfolding procedure. (see section 6-4)
3	Seat springs fall off	Please contact your dealer.
4	Seatback is not flipped down	Flip down the seatback and restart the folding / unfolding procedure.
5	Scooter is not on flat ground	Scooter cannot be folding/unfolding if it's not placed (lay down) on a flat ground. Please place the scooter on a flat ground and restart the folding / unfolding procedure.
6	Foreign Object on footplate.	Check and remove the object on the footplate and restart the folding / unfolding procedure.
7	Foreign Object underneath the seat	Check and remove the object underneath the seat and restart the folding / unfolding procedure.
8	Folding / unfolding light indicator on scooter is not in orange (auto version)	Auto version: please wake up the scooter. (see section 5)
	Folding / unfolding Release lever is stuck (manual version)	<ul style="list-style-type: none"> • Hold and pull the seat base backwards with one hand while • pulling the release lever with the other hand.
9	Other Error	The controller may have an internal fault. <ul style="list-style-type: none"> • Check all connections and wiring.

- If none of above correct problem, contact your authorized dealer. Do not disassemble the scooter by yourself.

9. OTHER INFORMATION

9.1 Options

We offer below options to add functionality for your scooter, please contact the authorized dealer for more information.

- We reserve the rights to modify or change the accessories.
- Please read through the cautions for battery and charger (section 6) if optional lead acid battery is purchased.
- Scooter folding / unfolding process may be interfered with optional accessories (armrests, bags) equipped. Please remove these accessories before folding/unfolding the scooter (accessories such as optional bigger battery will not affect the folding / unfolding, you don't have to remove the battery box).

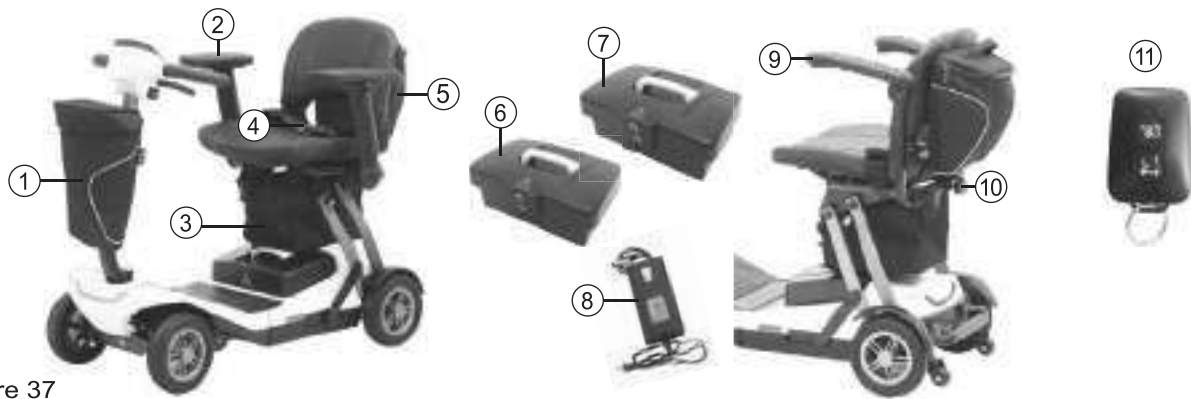

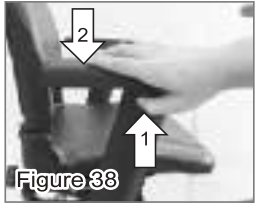



Figure 37

	Item	Qty	Description
01	Bag (Front)	1	Anti-splash, max. loading weight of 2Kgs. To be installed on the front tiller.
02	Adjustable Armrests	2	Multi-Adjustable with height and flip angles.
<div> <p>Angle adjustment</p> <ol style="list-style-type: none"> 1. Hold the armrest and pull up the latch. 2. Press downward the armrest. 3. Flip the armrest backward to the proper angle. <div>  <p>When flip the armrest backward, it is not in locking position. Please move the armrest back to the original position to lock.</p> </div> <div> <p>Height adjustment</p> <ol style="list-style-type: none"> 1. Press and hold the armrest latch. 2. Adjust the armrest height. 3. Release armrest latch to fix armrest in position. </div> <div>   </div> </div>			
03	Bag (underneath the seat)	1	Anti-splash, max. loading weight of 2Kgs. To be installed underneath the seat.
04	Seat Belt	1	To be installed for extra safety.
05	Bag (seatback)	1	Anti-splash, max. loading weight of 2Kgs. To be installed on the seatback.
06	Lead-acid Battery Pack	1	With lead-acid battery (12A x 2 pcs) comes with 24/2A lead-acid charger(Do NOT use it on lithium battery) Charging port for lead-acid charger : 3 ports
07	Lithium Battery Pack (bigger capacity 17.4AH)	1	With bigger 17.4Ah lithium battery, comes with 24V/4A lithium charger. (see item 8 for charger details)

08	Lithium Battery Charger (bigger capacity 24V/4A)	1	To be used on (17.4Ah x 1pcs) lithium battery pack, can be used on 11.6Ah lithium battery also. However, do NOT use it on lead-acid battery.
----	---	---	--

- Charging port for lithium charger : 4 ports
- Power Light : Red (power on)
- Charge Light : Red (normal charging) / Green(fully charged)



Figure 40

09	Compact Armrests	2	Compact armrests adjustable with flip angle.
10	Screws Knob for Armrests	2	For armrests fixation.
11	RF remote control (auto version):	1	A maximum of 4 remote controls can be set up on the same unit of scooter.

Remote Control External Frequency Pairing Setting

- Turn on the main power and unlock the scooter.
 1. Long-press the folding and unfolding buttons at the same time for 2 seconds. The light indicator (orange light) comes on to indicate it is ready for setup.
 2. Press either the folding or unfolding button once of the remote control.
 3. The light indicator will flash. When the indicator stops flashing, setting is complete.

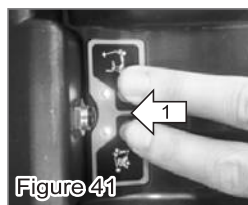


Figure 41



Change battery of remote control :

- Loosen the screw to open the back cover of remote control.
 1. Use screwdriver to loosen the two screws on the board.
 2. Take out the board and remove the battery (3V battery).
 3. Replace with new battery and follow the same steps backwards.


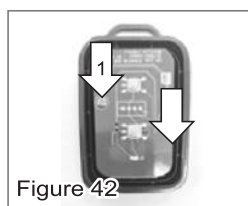


Figure 42

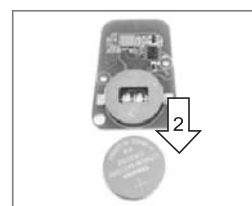

10. TECHNICAL SPECIFICATIONS

Specification	Li-Tech Velo
Overall Length	Folded : 775 mm / 30.5" Unfolded : 953 mm / 37.5"
Overall Width	Folded : 512 mm / 20.2" Unfolded : 512 mm / 20.2"
Overall Height	Folded : 540 mm / 21.3" Unfolded : 953 mm / 37.5"
Front Wheels	200 mm / 8"
Rear Wheels	200 mm / 8"
Weight W/ Batteries	31.2 kg / 68.9 lbs
Max. Speed	6.4 kmph / 4 mph
Weight Capacity	115 kg / 255 lbs
Ground Clearance	43 mm / 1.7"
Grade Climbable	6 degree
Curb Climbable	40mm / 1.6"
Turning Radius	1140 mm / 44.9"
Brake	Electro-Mechanical
Seat Type	Padded Foldable
Seat Width	419 mm / 16"
Motor Size	250W, 4900 r.p.m
Battery Size	24V. 17.4 Ah
Battery Weight	3.4 kg / 7.5 lbs
Travel Range	17 km / 10.6 Miles
Battery Charger	4A Off Board
Electronics	Sensor Key, Touch Screen, Coded Lock, Remote Control for Auto-Folding

*Subject to change without notice. (Rev. 2, 2023/01/13)

11. WARRANTY

Li-Tech warrants its Mobility Scooter product to be free from defects in materials and workmanship for a period of time from date of shipment, subject to certain limitations and exclusions as defined below.

Please read carefully the following terms and conditions to qualify for future repair service or replacement.

Frame: Li-Tech warrants the frame weldment of its Mobility Scooter to be free from defects in material and workmanship, under normal use for a period of 36 months for Mobility Scooters. Li-Tech will at its sole option, repair or replace any warrantable weldment or frame section, which fails within the warranty period.

Motor/Transaxle/Brake and Electronics: Li-Tech warrants the Motor/Transaxle/Brake/Electronics of its Mobility Scooter to be free from defects in material and workmanship, under normal use for a period of 24 months for Mobility Scooters. Li-Tech will at its sole option repair or replace any warrantable Motor/Transaxle/Brake or Electronic components, which fails within the warranty period.

WARRANTY EXCEPTIONS

- » Body panels, shrouds and cowlings
- » Consumable wear items such as: tyres, tubes, brake pads, suspension and motor brushes
- » Upholstery, seating and arm rests
- » Batteries carry a limited warranty of 7 months from the battery manufacturer
- » Any damage resulting from natural disaster
- » Any damage resulting from unauthorized repair or modification
- » Any damage resulting from negligent operation, storage or maintenance
- » Labor, service calls, shipping and any costs relating to loss of use

11. WARRANTY

NOTE

- All warrants replacement parts are warranted beginning on the shipment date of the mobility scooter or power chair when purchased.
- For a warranty claim, please submit the “Warranty Claim Form” with the information requested completed.
- Warranty is VOID if product is used in any sort of commercial or rental service.
- Warranty is VOID if any of the original serial numbers have been removed including from the controller, motor, transaxle, battery or battery charger.
- This warranty is for the replacement or repair, at the sole option of Li-Tech, for defective parts only. Li-Tech will cover the shipping charge, to our preference of freight, for only the replacement of part or repair of any such product and not other shipping charges that arise thereafter. Also note that this warranty does not cover any labor charges incurred during the replacement or repair.

THIS WARRANTY REPLACES AND SUPERCEDES ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE AND ANY OTHER OBLIGATIONS OR LIABILITIES ON THE PART OF LI-TECH, WHETHER IN CONTRACT, WARRANTY, NEGLIGENCE OR OTHERWISE. LI-TECH SHALL NOT BE LIABLE FOR AND HEREBY DISCLAIMS ALL CONSEQUENTIAL, INCIDENTAL AND CONTINGENT DAMAGES.

NOTES

 **LI-TECH**

velo

CareCo

1 Turing Court
Great Notley
Braintree
CM77 7AT